

Product Catalogue



incidental




MADE IN
TÜRKİYE

ENDOART®
FINAL TOUCH




ruby | dent



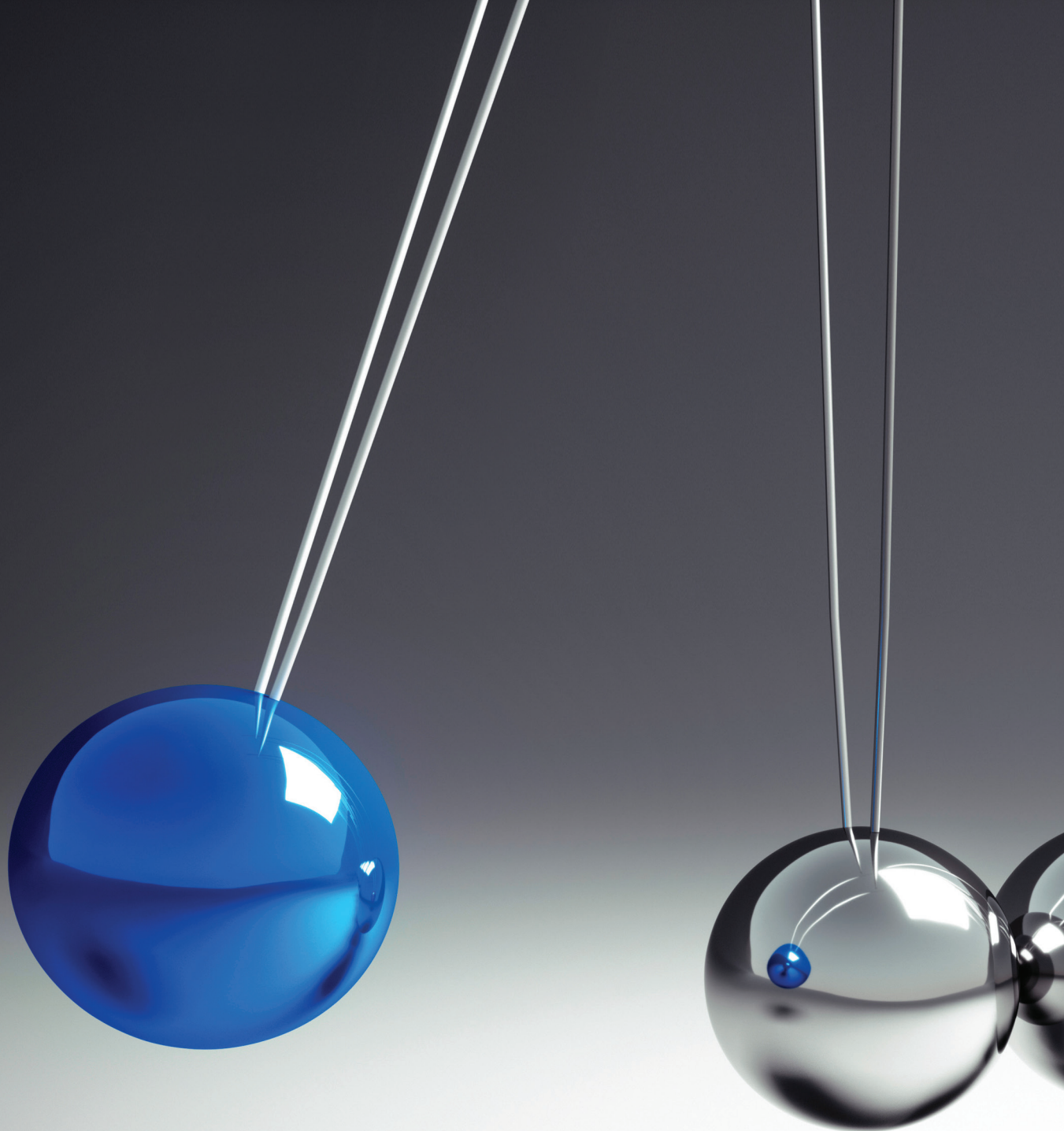

INDEX

INTRODUCTION	5	RubyLiner with Glass Ionomer	30
ENDODONTICS		RubyCal LC Calcium Hydroxide	30
EndoArt Action Rotary Files	7	Ruby Pit and Fissure Sealant	31
EndoArt Smart Rotary Files	8	FlexDam Gingiva Barrier	31
EndoArt Touch Rotary Files	9	RubyCem Luting	32
EndoArt Expert Reciprocating Files	11	RubyCem Filling	32
EndoArt Pedo Rotary Files	12	Ruby SE Cem	33
EndoArt Retreatment Files	13	Ruby Temp NE	33
EndoArt Path Rotary File	13	RubyZoe Zinc Oxide Eugenol Cement	34
EndoArt K Files	14	RubyZoe Eugenol	34
EndoArt H Files	14	RubyZoe Zinc Oxide Powder	34
EndoArt Spreader	14	RubyZinc P	35
EndoArt C Files	14	RubyZinc Poly	35
EndoArt Irrigation Tips	15	RubyShine Prophylaxis Paste	36
EndoArt Gutta Percha Points	16	Ruby Hemostat	37
EndoArt Absorbent Paper Points	16	Ruby Temp	37
EndoArt Endomotor Pro	17	RubyCap Capsule Amalgam	38
EndoArt Fiber Post	19	RubyCap NG Capsule Amalgam	39
EndoArt Bioceramic Sealer	20	RubyCap High Capsule Amalgam	39
EndoArt MTA Powder & Liquid	20	RubyGam Powder Amalgam	40
RubyCal Calcium Hydroxide Paste With Barium Sulphate	21	RubyGam High Powder Amalgam	40
RubyCal Calcium Hydroxide Paste With Iodoform	21	POLISHING & FINISHING	
RubyCal Calcium Hydroxide Powder & Liquid	21	RubyPlaton Polishing Discs	42
RubyCal Calcium Hydroxide Powder	21	RubyPlaton Polishing Strips	42
RubyEdta EDTA Gel	22	WEDGES	
EndoArt EDTA Solution	22	RubyPlaton Wooden Wedges	43
EndoArt Chlorex Gel	22	RubyPlaton Transparent Wedges	43
EndoArt Chlorex	22	ARTICULATING PAPERS	
EndoArt Sodium Hypochloride	23	RubyMark Articulating Paper Straight	44
EndoArt Citric Acid	23	RubyMark Articulating Paper Horseshoe	45
Ruby Iodoform Powder	23	RubyMark Articulating Paper 200µ Straight	45
Ruby Gutta Solver	23	AUXILIARY	
RESTORATIVES		Rubydent Curing Light	46
RubyFill Nano Composite	26	RubyCrownRemover	47
RubyComp Nano Composite	26	RubyPlaton Cotton Roll	47
RubyComp P Posterior Composite	27	ACRYLIC TEETH	
RubyFlow Flowable Composite	27	Compac Zahn	48
RubyFlow Bulk Fill Composite	27		
RubyBond Adhesive System	28		
RubySEBond Self Etching Adhesive System	29		
RubyEtch Etching Gel	29		
RubyEtch HF Porcelain Etch	29		

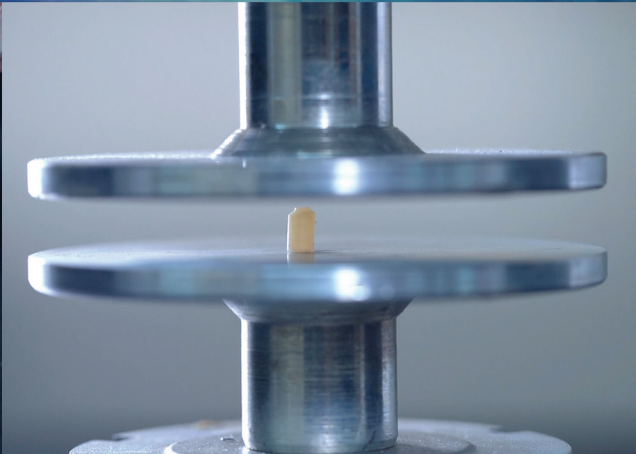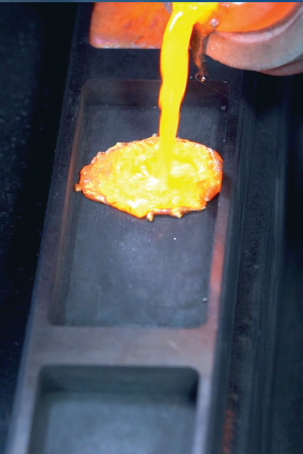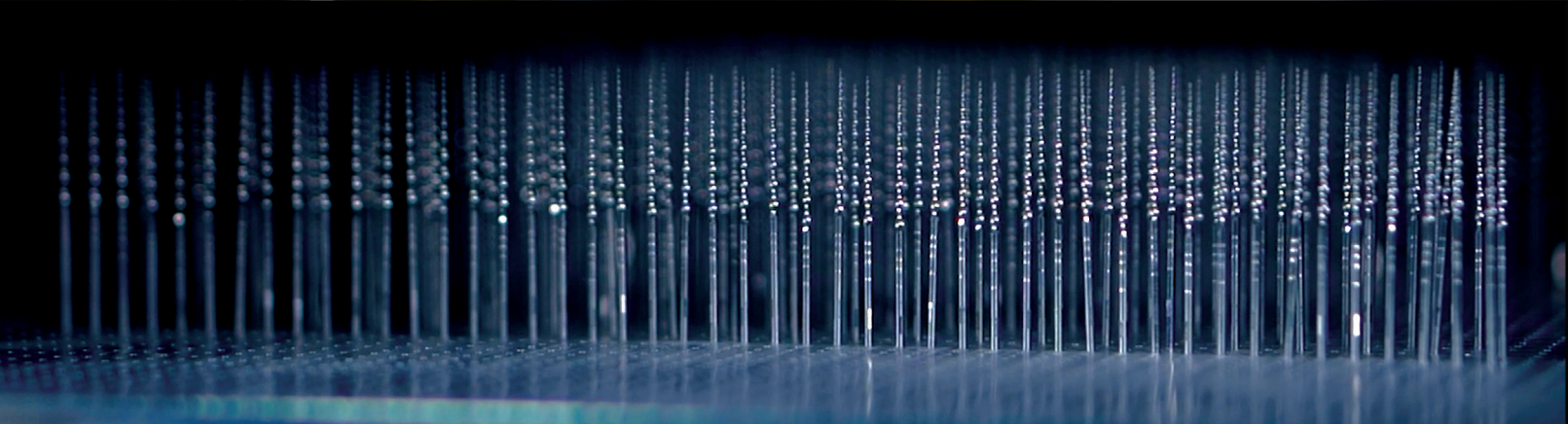

About us

Company

Inci Dental was founded 1994 in Istanbul. Since then the company produce a variety of consumable materials. Since 1994, the concept of quality is established, always giving priority to customer satisfaction, the production process by making significant investments in R & D activities, cooperate with universities and Scientific and Technological Research Council of Turkey, made production as newest and unique products, in the sector so has achieved a leading place.

The company is only producer of variety materials in Türkiye. Inci Dental is manufacturer of Dental amalgam filling material, dental composites filling material, endodontic rotary files, fiber posts, polishing discs, articulating paper, acrylic teeth and etc. We export our products more than 80 countries in the World.

Our quality assurance approve by European Union with the certificates of ISO 9001, ISO 13485 and the EC.

Vision

To be a company which is the first choice of domestic and foreign business partners.

In this context, the sale of the products that we have done all, are the first brands that come to mind of the dentists.

Mission

Inci dental, in dental sector;

Who makes the difference

Who find, research and use the right technologies in all over the world

Who adapt quickly to technological developments

Who produce best quality

Who lead the innovation

A company which holds the highest level of employee, domestic and foreign business partners satisfaction.


















incidental















ENDOART®
FINAL TOUCH



**Heat Treated Control Memory
Rotary Files**

SX	S1	S2	F1	F2	F3	F4	F5
							
							
Starter	Shaping						
Torque: 3.0 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm



SX	S1	S2	F1	F2	F3	F4	F5
							
							
Starter	Shaping						
Torque: 2.6 Ncm 300 Rpm	Torque: 1.8 Ncm 300 Rpm	Torque: 1.8 Ncm 300 Rpm	Torque: 1.8 Ncm 300 Rpm	Torque: 1.8 Ncm 300 Rpm	Torque: 1.8 Ncm 300 Rpm	Torque: 1.8 Ncm 300 Rpm	Torque: 2.6 Ncm 300 Rpm

EndoArt Action rotary files allow you to move comfortably in the canal thanks to the special structure of its surface, while its variable taper is designed to evacuate debris from the canal. The developed special heat treatment technology ensures that the Ni-Ti wire has a much more flexible and durable structure and also turns it into a Gold or Blue color. EndoArt Action Gold files have a higher sharpness efficiency and higher fracture resistance than standard Ni-Ti files. EndoArt Action Blue files allows pre-bending before entering the canal thanks to its high flexibility efficiency with special control memory and it has a much higher fracture resistance than standard Ni-Ti files. It adapts to the shape of the canal in curvature structures.

Action Kit SX - S1 - S2 - F1 - F2 - F3

Action Gold Rotary Files

Rotation Files
Variable Taper
Better Cutting Efficiency
x2 times better Cyclic Fatigue Resistance than conventional Files
Excellent Strength
No bounce back
Respect Original Anatomy of Root Canal
Great Fracture Resistance
Lasered Taper and Size information on the shank
Autoclavable Box

Control Memory
Heat Treated Wire

Convex Triangle
Cross - Section



Action Blue Rotary Files

Rotation Files
Variable Taper
Better Flexibility
X4 times better Cyclic Fatigue Resistance than conventional Files
Excellent Strength
No bounce back
Respect Original Anatomy of Root Canal
Great Fracture Resistance
Lasered Taper and Size information on the shank
Autoclavable Box

Wire : Gold
Length : 21mm, 25mm, 31mm
Packaging : 6 Files



Wire : Blue
Length : 21mm, 25mm, 31mm
Packaging : 6 Files



WhatsApp
Get information!

SMART



10/04	15/04	20/04	25/04	30/04	35/04	40/04	50/04
04 Taper		Starter	Shaping / Finishing	Finishing			
Torque: 3.0 Ncm 300 Rpm		Torque: 2.0 Ncm 300 Rpm	Torque: 2.0 Ncm 300 Rpm	Torque: 2.0 Ncm 300 Rpm	Torque: 2.0 Ncm 300 Rpm	Torque: 2.0 Ncm 300 Rpm	Torque: 2.0 Ncm 300 Rpm



10/04	15/04	20/04	25/04	30/04	35/04	40/04	50/04
04 Taper		Starter	Shaping / Finishing	Finishing			
Torque: 3.0 Ncm 300 Rpm		Torque: 1.6 Ncm 300 Rpm	Torque: 1.6 Ncm 300 Rpm	Torque: 1.6 Ncm 300 Rpm	Torque: 1.6 Ncm 300 Rpm	Torque: 1.6 Ncm 300 Rpm	Torque: 1.6 Ncm 300 Rpm

EndoArt Smart rotary files have fixed %4 taper allow you to work more safely in narrow canals. The developed special heat treatment technology ensures that the Ni-Ti wire has a much more flexible and durable structure and also turns it into a Gold or Blue color. In addition to its high sharpness efficiency, EndoArt Smart Gold has a higher fracture resistance than standard Ni-Ti files. EndoArt Smart Blue allows pre-bending before entering the canal thanks to its high flexibility efficiency with special control memory and it has a much higher fracture resistance than standard Ni-Ti files. It adapts to the shape of the canal in cases with high slopes.

Smart Kit 15/04-20/04-25/04-30/04-35/04-40/04 19/08-10/04-15/04-20/04-25/04-30/04 19/08-15/04-20/04-25/04-30/04-35/04 19/08-16/02-20/04-25/04-30/04-35/04

Smart Gold Rotary Files

- Rotation Files
- Regular Taper
- Better Cutting Efficiency
- x2 times better Cyclic Fatigue Resistance than conventional Files
- Excellent Strength
- No bounce back
- Respect Original Anatomy of Root Canal
- Great Fracture Resistance
- Lasered Taper and Size information on the shank
- Autoclavable Box



Smart Blue Rotary Files

- Rotation Files
- Regular Taper
- Better Flexibility
- X4 times better Cyclic Fatigue Resistance than conventional Files
- Excellent Strength
- No bounce back
- Respect Original Anatomy of Root Canal
- Great Fracture Resistance
- Lasered Taper and Size information on the shank
- Autoclavable Box


Wire : Gold
Length : 21mm, 25mm, 31mm
Packaging : 6 Files



Wire : Blue
Length : 21mm, 25mm, 31mm
Packaging : 6 Files



WhatsApp
Get information!

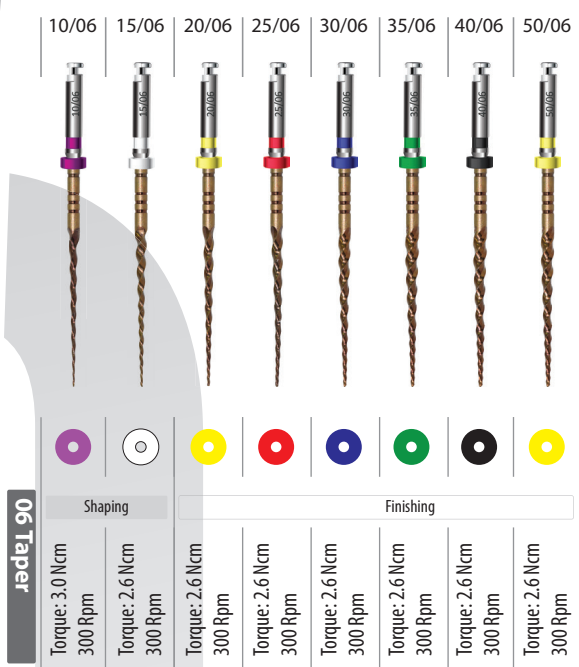

EndoArt Touch rotary files fixed %6 taper allow you to prepare at a standard angle along the canal. The developed special heat treatment technology ensures that the Ni-Ti wire has a much more flexible and durable structure and also turns it into a Gold or Blue color. In addition to its high sharpness efficiency, EndoArt Touch Gold has a higher fracture resistance than standard Ni-Ti files. EndoArt Touch Blue allows pre-bending before entering the canal thanks to its high flexibility efficiency with special control memory and compared to standard Ni-Ti files, it has a much higher fracture resistance. It adapts to the shape of the canal in cases with high slopes.

Touch Kit 15/06-20/06-25/06-30/06-35/06-40/06 19/08-15/04-20/04-20/06-25/06-30/06 19/08-10/06-15/06-20/06-25/06-30/06 19/08-15/06-20/06-25/06-30/06-35/06

Touch Gold Rotary Files

Rotation Files
Regular Taper
Better Cutting Efficiency
x2 times better Cyclic Fatigue Resistance than conventional Files
Excellent Strength
No bounce back
Respect Original Anatomy of Root Canal
Great Fracture Resistance
Lasered Taper and Size information on the shank
Autoclavable Box

Control Memory
Heat Treated Wire

Convex Triangle
Cross - Section



Touch Blue Rotary Files

Rotation Files
Regular Taper
Better Flexibility
X4 times better Cyclic Fatigue Resistance than conventional Files
Excellent Strength
No bounce back
Respect Original Anatomy of Root Canal
Great Fracture Resistance
Lasered Taper and Size information on the shank
Autoclavable Box

Wire : Gold
Length : 21mm, 25mm, 31mm
Packaging : 6 Files



Wire : Blue
Length : 21mm, 25mm, 31mm
Packaging : 6 Files



WhatsApp
Get information!

ENDART[®]

FINAL TOUCH



Reciprocation




EndoArt Expert files work with a reciprocating motion that allows you to perform root canal treatment with a single file. It includes alternative tapers, %6 and %4 angles for single-file systems, and allows you to prepare at a standard angle along the canal. EndoArt Expert files, which can be used with endomotors capable of reciprocating movement, work with a rotation angle adjustment of 150° left and 30° right. The developed special heat treatment technology makes the Ni-Ti wire a much more flexible and durable structure. EndoArt Expert Gold files have higher sharpness and higher fracture resistance than standard Ni-Ti files. EndoArt Expert Blue allows pre-bending before entering the canal thanks to its high flexibility efficiency with special control memory and compared to standard Ni-Ti files, it has a much higher fracture resistance.

Expert Kit R25-R25-R25-R25-R40-R50 R25-R25-R40-R40-R50-R50 R25-R25-R40-R50 R25-R40-R50

Expert Gold Reciprocation Files

Reciprocation Files
Regular Taper
Better Cutting Efficiency
x2 times better Cyclic Fatigue Resistance than conventional Files
Excellent Strength
No bounce back
Respect Original Anatomy of Root Canal
Great Fracture Resistance
Lasered Taper and Size information on the shank
Autoclavable Box

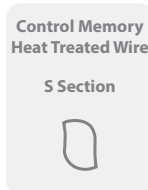

Wire : Gold
Length : 21mm, 25mm, 31mm
Packaging : 6 Files

Expert Blue Reciprocation Files


Reciprocation Files
Regular Taper
Better Flexibility
X4 times better Cyclic Fatigue Resistance than conventional Files
Excellent Strength
No bounce back
Respect Original Anatomy of Root Canal
Great Fracture Resistance
Lasered Taper and Size information on the shank
Autoclavable Box

Wire : Blue
Length : 21mm, 25mm, 31mm
Packaging : 6 Files





PEDO

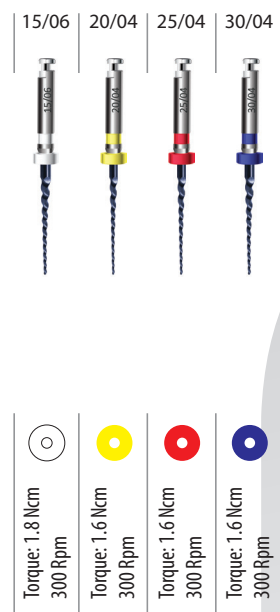

EndoArt Pedo rotary files are produced specially for pediatric cases and has a length of 18 mm. The improved heat treatment technology of EndoArt Pedo ensures that the Ni-Ti wire has a much more flexible and durable structure. EndoArt Pedo files have higher sharpness and fracture resistance than standard Ni-Ti files.

Pedo Kit 15/06-20/04-25/04 15/06-20/04-30/04

Pedo Gold Rotary Files

- Rotation Files
- Specially produced for pedodontic cases
- Regular Taper
- Better Cutting Efficiency
- x2 times better Cyclic Fatigue Resistance than conventional Files
- Excellent Strength
- No bounce back
- Respect Original Anatomy of Root Canal
- Great Fracture Resistance
- Lasered Taper and Size information on the shank
- Autoclavable Box



Wire	: Gold
Length	: 18mm
Packaging	: 3 Files

Pedo Blue Rotary Files

- Rotation Files
- Specially produced for pedodontic cases
- Regular Taper
- Better Flexibility
- X4 times better Cyclic Fatigue Resistance than conventional Files
- Excellent Strength
- No bounce back
- Respect Original Anatomy of Root Canal
- Great Fracture Resistance
- Lasered Taper and Size information on the shank
- Autoclavable Box

Wire	: Blue
Length	: 18mm
Packaging	: 3 Files



WhatsApp
Get information!

RETREATMENT



EndoArt Retreatment rotary files developed for repeated root canal treatments. Thanks to its high cutting efficiency, it is used in the process of removing the materials in the canals. The system consists of three instruments, D1 - D2 and D3 which are designed for cleaning the canal filling in the (D1) coronal, (D2) middle, and (D3) apical thirds. While D1 has an active tip at the first entry into the root canal, the tips of continuing instruments are not cutting.

Retreatment Kit D1-D1-D2-D2-D3-D3

Retreatment Files

- Rotation Files
- Specially produced for retreatment cases
- Designed to remove filling materials like Gutta Percha
- x2 times better Cyclic Fatigue Resistance than conventional Files
- Better Cutting Efficiency
- Great Fracture Resistance
- Lasered Taper and Size information on the shank
- Autoclavable Box

Wire	: Gold
Length	: 16mm, 18mm, 22mm
Packaging	: 6 Files

Control Memory
Heat Treated Wire

Convex Triangle
Cross - Section

PATH FILE



Path files are used to create a guide (Glide) path. The use of angled files in combination with hand files before use, reduces the risk of step formation and perforation. The size and taper of the file are 16/02. The ends of the file are rounded to prevent the formation of step and zip formation.

Path Gold Rotary Files

- Rotation Files
- Specially produced for creating glide path
- Easily reach apex
- Better Flexibility
- x2 times better Cyclic Fatigue Resistance than conventional Files
- Great Fracture Resistance
- Lasered Taper and Size information on the shank
- Autoclavable Box

Wire	: Gold
Length	: 21mm, 25mm, 31mm
Packaging	: 6 Files

Control Memory
Heat Treated Wire

Quadrangular





EndoArt K Files

Codes: 06 08 10 15 20 25 30 35 40 45
50 55 60 70 80 90 100 110 120 130 140

Ass. 15/40 – Ass. 45/80 – Ass. 90/140

Stainless Steel

Length: 21 mm, 25 mm, 31 mm

6 pieces in a box



EndoArt H Files

Codes: 06 08 10 15 20 25 30 35 40 45
50 55 60 70 80 90 100 110 120 130 140

Ass. 15/40 – Ass. 45/80 – Ass. 90/140

Stainless Steel

Length: 21 mm, 25 mm, 31 mm

6 pieces in a box



EndoArt Spreader

Codes: 15 20 25 30 35 40/ Ass. 15/40

Stainless Steel

Length: 21 mm, 25 mm

6 pieces in a box



EndoArt C Files

Codes: 06 08 10 12.5 15/ Ass. 06-10

Stainless Steel

Length: 25 mm

6 pieces in a blister package



Get information!



EndoArt Irrigation Tips

It is plastic tip for irrigation during root canal treatment. It provides effective cleaning and disinfection thanks to double side vent. It helps to remove debris and smear layer. It has conical shape to adapt root canal anatomy. Made by special polypropylen material.

Properties

30G – 0.3 mm tip
Can be autoclavable

Packaging

20 pcs in a box



WhatsApp
Get information!

EndoArt Gutta Percha Points

Gutta Percha root canal filling material for easy fillings. It has variable taper and sizes.

Properties

02/04/06 Taper Options
Excellent Radiopacity
Length Marked
Hand Rolled



EndoArt Paper Points

Absorbent Paper Point to dry root canals.

Properties

Highly Absorbent
Sterile Packaging
Variety Sizes and Tapers
Length Marked
Hand Rolled





EndoArt Endomotor Pro

Cordless endomotor with apex locator function.

Properties

Ready program settings for different files
 3 Rotation Mode; Clockwise/Anti Clockwise and Reciprocating
 Min 50 Rpm - Max 1800 Rpm and 5.0 N.Cm
 Micro USB Charging
 16:1 Contra Angle
 Auto Start in the root canal
 Auto Reverse
 Brushless Motor
 360 ° rotate centre angle
 Plastic centre angle option



6:1

OLED
screen



AUTO REVERSE



AUTO START-STOP



Reciprocation



WhatsApp
Get information!

ENDOART®

FINAL TOUCH

Fiber Post





EndoArt Fiber Post

Universal resin glass fiber post system for esthetic and metal free restorations. It has different type of diameters which are 1,2 mm, 1,4 mm, 1,6 mm and 1,8 mm.

Properties

Radiopaque
4 drills available inside of the kit
Conical tip
High resistance
Colour coded
It is compatible with adhesive resin cements

Packaging


EndoArt Fiber Post Kit – 5 pcs. 1,2 mm posts, 5 pcs. 1,4 mm posts, 5 pcs. 1,6 mm posts, 5 pcs. 1,8 mm posts, 4 drills assorted
Ø 1.2 - Ø 1.4 - Ø 1.6 - Ø 1.8

EndoArt Fiber Post Refill - 5 pcs. 1,2 mm posts

EndoArt Fiber Post Refill - 5 pcs. 1,4 mm posts

EndoArt Fiber Post Refill - 5 pcs. 1,6 mm posts

EndoArt Fiber Post Refill - 5 pcs. 1,8 mm posts

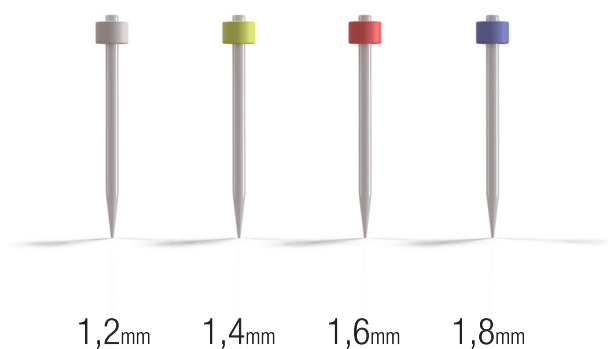

1,2mm

1,4mm

1,6mm

1,8mm





EndoArt Bioceramic Sealer

It is highly biocompatible premixed injectable calcium aluminosilicate paste for root canal filling and sealing. It provides superior root and obtention and sealing.

Properties

Resin free	Highly flowable
High pH	Non toxic
Zero shrinkage	Good Setting Time
Radiopaque	No mixing

Packaging

2 gr in a syringe



EndoArt MTA Powder&Liquid

MTA (Mineral Trioxide Aggregate) is tricalcium silicate-based root canal repair material. Biocompatible and fine grain powder allows you to make easy and successful treatment. It can be used for root canal perforation, apical plug during apexification, pulp capping, pulpotomy and root canal filling material.

Properties

Excellent Handling Properties
High pH
Promotes mineralization of tissue
Radiopaque
No discolouration

Packaging

0,14 gr MTA powder - 5 ml Liquid
0,14 gr MTA powder x 7 - 5 ml Liquid
1 gr MTA Powder – 5 ml Liquid



RubyCal Calcium Hydroxide Paste With Barium Sulphate

Calcium hydroxide premixed paste non-setting for root canal therapy and vital pulp therapy. It comes with application tips.

Properties

High Radiopacity
Easy removal
High pH > 12
Excellent handling characteristics

Packaging

2 gr syringe
2 gr syringe x 2
2 gr syringe x 4



RubyCal Calcium Hydroxide Paste With Iodoform

Calcium hydroxide premixed paste non-setting for root canal therapy and vital pulp therapy. It comes with application tips.

Properties

Radiopacity
Promote apexification
Disinfects Canal
Excellent handling characteristics

Packaging

2 gr syringe
2 gr syringe x 2
2 gr syringe x 4



RubyCal Calcium Hydroxide Powder & Liquid

Calcium hydroxide powder and liquid for root canal therapy

Properties

Temporary Root Filling
High pH > 12
Excellent handling characteristics
Direct and indirect pulp capping
Base under non-adhesive restorations
Easy to mix

Packaging

10 gr powder and 10 ml liquid



RubyCal Calcium Hydroxide Powder

Calcium hydroxide powder for root canal therapy.

Properties

Temporary Root Filling
Direct and indirect pulp capping
Base under non-adhesive restorations
Antibacterial Properties
Easy to mix

Packaging

30 gr powder



RubyEdta 17% EDTA Gel

It is 17% EDTA based root canal preparation gel for lubricating and cleaning of root canal. It comes with application tips.

Properties

Easily moving
Antibacterial effect
Reduce file breakage

Packaging

3 ml syringe
3 ml syringe x 2
3 ml syringe x 4



EndoArt Chlorex Gel

EndoArt Chlorex is %2 chlorhexidine gel for preparation of root canal.

Properties

2% chlorhexidine digluconate
Easier root canal treatment
Minimum risk of root canal breakage

Packaging

2 gr in a syringe



EndoArt EDTA Solution

EndoArt EDTA is irrigation solution for preparation and cleaning of root canals, remove smear layer. %17 EDTA or %5 EDTA options are available.

Properties

Clean root canals rapidly
Neutral PH
Helps to open calcified canals

Packaging

250 ml liquid
100 ml liquid



EndoArt Chlorex

EndoArt Chlorex is %2 chlorhexidine digluconate irrigation liquid for root canal treatment.

Properties

2% chlorhexidine digluconate
Improves mechanical preparation and cleaning of root canals

Packaging

250 ml liquid
100 ml liquid



EndoArt Sodium Hypochlorite

It is 5.25% sodium hypochlorite irrigation liquid for root canal treatment.

Properties

It has surfactant increase penetration to dentinal tubules.
Removes smear layer
Effective canal cleaning
Strong antibacterial effect

Packaging

250 ml liquid



EndoArt Citric Acid

It is %40 citric acid irrigation liquid for root canal treatment.

Properties

Removes smear layer from the dentinal walls
Effective canal cleaning
Cleaning calcium hydroxides from the root canals
Improves the mechanical cleaning of the root canal

Packaging

250 ml liquid



Ruby Iodoform Powder

Iodoform is used for root canal treatment. Can be used for filling of root canals after depulping, chronic periodontitis, wound dressing.

Properties

Pure Iodoform (Pharma Grade)
High PH 12,5

Packaging

30 gr powder



Ruby Gutta Solver

RubyGuttaSolver is used for removing of gutta percha materials during re-treatment of root canals.

Properties

%100 Orange Oil
Non toxic, natural
Easy to wash out

Packaging

10 ml liquid

ruby | dent

ENDART®
FINAL TOUCH



MADE IN
TÜRKİYE



incidental



RubyFill Nano Composite

RubyFill Nano is radiopaque nano hybrid composite system which can be applied in the anterior and posterior area.

Properties

Universal Nano Composite
It has light refraction property like a normal tooth due to the nano fillers in it
It has high filler rate
It can be polished high with high degree
It has high resistance to abrasion
It provides ease of use

Packaging

Refill – 4 gr in a syringe
Kit 1 - 4 gr syringe x 6 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
Shades – A1 – A2 – A3 – A3,5 – B1 – B2 – C1 – D3 – OA2 – Extra Bleach



RubyComp Nano Composite

RubyComp Nano is radiopaque nano hybrid composite system which can be applied in the anterior and posterior field.

Properties

Higher Esthetic
Low shrinkage
It has high filler rate
It can be polished with high degree
It has high resistance to abrasion
It provides ease of use
Long durability
Non sticky

Packaging

Refill – 4 gr in a syringe
Kit 1 - 4 gr syringe x 6 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
Shades – A1 – A2 – A3 – A3,5 – B1 – B2 – C1 – D3 – OA2 – Extra Bleach



ISO 4049

CLASS 2A



 **WhatsApp**
Get information!



RubyComp P Posterior Composite

Rubycomp P is composite system which can be applied in the posterior field.

Properties

Non stick to instruments
It has surface hardness is like real enamel
It has high intensity
It has a perfect consistence which can be shaped/formed easily
It has high filler rate
It resists to chewing force
It does not discolor after the polimeration.

Packaging

Refill – 4 gr in a syringe
Kit 1 - 4 gr syringe x 6 pcs. , 5 ml RubyBond x 1 pcs. 3 ml RubyEtch x 1 pcs. and Accessories
Shades – A1 – A2 – A3 – A3,5



RubyFlow Flowable Composite

RubyFlow is flowable light curing radiopaque composite.

Properties

Low shrinkage, low wear
High gloss polishability
Beatiful esthetic restorations
Very low surface roughness
Excellent handling
Non sticky
Very good marginal adaptation

Packaging

Refill – 2 gr in a syringe
Kit 1 - 2 gr in a syringe x 2 pcs.
Kit 2 - 2 gr syringe x 4 pcs.
Shades – A1 – A2 – A3 – A3,5 – B1 - OA2

RubyFlow Bulk Fill

Flowable bulk-fill composite for class I and II restorations can be applied on the posterior field using bulk fill technique. Thanks to material flowability the material adapts perfectly to the cavity walls.

Properties

Filling up to 4 mm
Different colour options
Low shrinkage
Increased mechanical properties
Single Layer
High Radiopacity

Packaging

2 gr in a syringe



ISO 4049

CLASS 2A



 **WhatsApp**
Get information!

5th.
Generation



RubyBond Adhesive System

RubyBond is light curing single component 5th generation bonding agent for the adhesive restorations. It is universal and compatible with all composites.

Properties

High bond strength
Perfect adhesion to dentin and enamel
5th generation

Packaging

5 ml in a bottle
5 ml in a bottle + 3ml Etching

7th.
Generation



RubySEBond Self Etching Adhesive System

RubySEBond is self-etching light curing one component 7th generation adhesive for the adhesive restorations.

Properties

Very strong bonding to enamel and dentine
No extra etching with phosphoric acid
Long lasting restorations
Time savings
7th generation

Packaging

5 ml in a bottle



CLASS 2A



WhatsApp
Get information!



RubyEtch Etching Gel

RubyEtch is an etching gel containing 37% phosphoric acid. Delivered in 3 ml. and 50 ml syringe with very thin application needles.

Properties

Easily Rinsed
Blue coloured
Optimum Surface Etching

Packaging

3 ml syringe
3 ml syringe x 2
3 ml syringe x 4
50 ml syringe



RubyEtch HF Porcelain Etch

RubyEtch HF is 9.5% Buffered Hydrofluoric acid in gel form. Used to improve the mechanical bonding between porcelain and adhesive agent.

Properties

Porcelain Restoration Repair
Buffered Formula
Easily Rinsed

Packaging

3 ml syringe
Porcelain Repair Kit - 3 ml syringe + Silane



CLASS 2A



 **WhatsApp**
Get information!



RubyLiner with Glass Ionomer

Light curing glass ionomer based cavity liner to protect pulp of the tooth. It is used under direct and indirect restorations, like composites, amalgams, ceramics and another metal restorations.

Properties

- Excellent biocompatibility
- Releases fluoride
- Radiopaque
- Excellent handling characteristics

Packaging

- 2 gr syringe
- 2 gr syringe x 2
- 2 gr syringe x 4



RubyCal LC Calcium Hydroxide

Light curing calcium hydroxide for pulp capping and liner material under dental filling.

Properties

- Direct & indirect pulp capping
- Light curing
- Single component

Packaging

- 2 gr syringe
- 2 gr syringe x 2
- 2 gr syringe x 4



ISO 4049

CLASS 2A



[WhatsApp](#)
Get information!



Ruby Pit and Fissure Sealant

Light curing fluoride releasing sealant for occlusal surfaces. It can be used as preventive products.

Properties

Fluoride releasing
Radiopaque
Easy application
Low shrinkage

Packaging

2 gr syringe
2 gr syringe x 2
2 gr syringe x 4



FlexDam Gingiva Barrier

Light curing resin-based gingiva barrier to protect soft tissue.

Properties

Easily removable after treatment
It minimizes heat and soft tissue sensitivity during curing
Blue shade
5 second curing

Packaging

2 gr syringe
2 gr syringe x 2
2 gr syringe x 4



CLASS 2A



WhatsApp
Get information!



RubyCem Luting

Biocompatible radiopaque powder and liquid glass ionomer luting cement. It is used for cementation of inlays, onlays, crowns and bridges made of metal or porcelain fused to metal.

Properties

- Resin free
- High flexural strength
- High fluoride release
- Radiopaque
- Easy mixing

Packaging

10 gr powder and 10ml liquid



RubyCem Filling

Biocompatible radiopaque powder and liquid glass ionomer filling cement.

Properties

- Resin free
- Great aesthetics
- Radiopaque
- Easy mixing
- Very Durable
- Higher compressive strength

Packaging

10 gr powder and 10ml liquid



Ruby SE Cem

Self-etching and self-adhesive permanent resin cement for luting. It is used for cementation of: crowns, bridges, inlays, onlays and posts.

Properties

- Resin based
- Higher bonding strength
- Excellent marginal sealing
- Time saving
- High performance



Ruby Temp NE

Zinc oxide non eugenol temporary cement. It comes with automix syringe and mixing tips. It is used for temporary cementation of crowns, bridges. It can be used for patients showing allergic reaction to the eugenol.

Properties

- Free of eugenol
- Easy to use
- Perfect sealing
- Low solubility



RubyZoe Zinc Oxide Eugenol Cement

Zinc oxide and eugenol based dental cement for temporary filling.

Properties

Temporary Filling
Root canal sealing
Root canal filling

Radiopaque
Base and lining material under fillings
Easy to mix

Packaging

45gr powder and 15 ml liquid



RubyZoe Eugenol

Eugenol suitable for use with Zinc oxide powder for temporary filling.

Properties

%100 Eugenol (Pharma Grade)
Easy to mix

Packaging

15 ml liquid



RubyZoe Zinc Oxide Powder

Zinc oxide powder mixing with eugenol for temporary filling.

Properties

%100 Zinc Oxide (Pharma Grade)
Base and lining material under fillings
Temporary Fixing crown and bridges
Indirect pulp capping
Easy to mix

Packaging

45gr powder



Ruby Zinc P

Zinc phosphate cement for crowns and bridges and as liner material for all type of filling material.

Properties

- High compressive strength
- Radiopaque
- Easy application
- Low film thickness



Ruby Zinc Poly

Zinc polycarboxylate cement for cementation of Crowns, bridges, inlays, onlays. It is also used as Base/lining under composite, amalgam or glass-ionomer and as temporary filling material.

Properties

- High adhesive strength
- Contains no phosphoric acid
- Radiopaque
- Easy mixing and application
- Low film thickness



RubyShine Prophylaxis Paste

RubyShine Prophylaxis Paste is specially formulated of polishing and cleaning agents which has nano particuls in it. Available in mint, strawberry, bubblegum and green apple flavor.

Packaging

30 gr in a bottle
150 gr in a bottle
340 gr in a bottle
2 gr x 100 pcs. in a box
2 gr x 200 pcs. in a box



Strawberry



Mint



Bubblegum



Orange



Class I



 **WhatsApp**
Get information!



RubyHemostat

25% buffered aluminium sulphate solution. It is used for stopping minor gingival bleeding. It can be used during restorative procedures, crown preparations and placement of bands, matrixes and before impression taking.

Properties

Epinephrine free
Stable
Quick hemostasis

Packaging

10ml liquid



RubyTemp

Temporary filling material for proven filling and temporary insertion of permanent crowns and bridges

Properties

Good marginal adaptation
Self curing with action of saliva
No mixing
Non eugenol

Packaging

30 gr in bottle



RubyCap Capsule Amalgam

High copper alloy which contains a mixture of spherical and irregular particles. Predosed in self activating capsules. Gamma-2 free. It has 4 spill options as spill 1, spill 2, spill 3 and spill 5.

Indications

Stress bearing fillings in posterior teeth (Class I and II)
Retrograde root canal filling after apex resection

Features

Mixture of spherical and irregular particles
High copper alloy
High affinity with mercury
Firm condensation
Smooth carving
Excellent condensability

Composition

45% Ag – 30,5% Sn – 24% Cu – 0,5% Zn
45% Ag – 31% Sn – 24% Cu (Zinc Free)

Packaging

50 capsule in a jar



RubyCap NG Capsule Amalgam

High copper alloy which contains a mixture of spherical and irregular particles. It contains higher silver (55% Ag) than RubyCap. Predosed in self activating capsules. Gamma-2 free. It has 4 spill options as spill 1, spill 2, spill 3 and spill 5.

Indications

Stress bearing fillings in posterior teeth (Class I and II)
Retrograde root canal filling after apex resection

Features

High silver / high copper alloy
Easy to use

Composition

55% Ag – 29% Sn – 15,8% Cu – 0,2% Zn

Packaging

50 capsule in a jar



ISO 20749

CLASS 2A



WhatsApp
Get information!



RubyCap High Capsule Amalgam

High silver alloy which contains a mixture of spherical and irregular particles. Predosed in self activating capsules. Gamma-2 free. It has 4 spill options as spill 1, spill 2, spill 3 and spill 5.

Indications

Stress bearing fillings in posterior teeth (Class I and II)
Retrograde root canal filling after apex resection

Features

Mixture of spherical and irregular particles
High silver content
Exceptionally durable
Easy to process
Relatively fast hardening

Composition

%69,2 Ag %18,6 Sn %11,9 Cu %0,3 Zn
%69,2 Ag %18,9 Sn %11,9 Cu (Zinc Free)

Packaging

50 capsule in a jar

Technical Properties

A well balanced alloy for quality fillings;
Spill 1 – 400 mg Alloy / 360 mg Mercury
Spill 2 – 600 mg Alloy / 540 mg Mercury
Spill 3 – 800 mg Alloy / 720 mg Mercury
Spill 5 – 1200 mg Alloy / 1080 mg Mercury

Working Time: 4 minutes
Carving & burnishing: 5 – 9 minutes

Conforms to:
ISO 20749 Pre-capsulated Dental Amalgam

Compressive Strength;
1h < 397 MPa
24h < 437 MPa

Creep : %0,99
Dimensional Change During Hardening : %0,01





RubyGam Powder Amalgam

High copper alloy which contains a mixture of spherical and irregular particles. Gamma-2 free.

Indications

Stress-bearing fillings in posterior teeth (Class I and II)
Retrograde root canal filling after apex resection

Features

Mixture of spherical and irregular particles
High copper alloy
High affinity with mercury
Firm condensation
Smooth carving
Excellent condensability

Composition

45% Ag – 30,5% Sn – 24% Cu – 0,5% Zn

Packaging

Boxes of 30 gr



RubyGam High Powder Amalgam

High silver alloy which contains a mixture of spherical and irregular particles. Gamma-2 free.

Indications

Stress bearing fillings in posterior teeth (Class I and II)
Retrograde root canal filling after apex resection

Features

A mixture of spherical and irregular particles
High silver content
Exceptionally durable
Easy to process
Relatively fast hardening

Composition

69,2% Ag – 18,6% Sn – 11,9% Cu – 0,3% Zn

Packaging

Boxes of 30 gr



ISO 20749

CLASS 2A



WhatsApp

Get information!

ruby | dent



Ruby Cap

STABILITY

HIGH STRENGTH

EXCELLENT HANDLING

NON GAMMA 2

RubyPlaton Polishing Discs

Has different kind of discs which has different grins and different sizes. It allows excellent results in polishing. Flexible and easy to use

Properties

14 mm blue, yellow, white, green, brown, pink

- **White** – Extra Fine for Polishing
- **Yellow** – Fine for Finishing
- **Green** – Medium for Countoring
- **Blue** – Coarse for Gross Reduction
- **Brown** – X Coarse for Gross Reduction
- **Pink** – XX Coarse for Gross Reduction

Packaging

10 pcs discs and 1 mandren in a pack
 40 pcs discs and 1 mandren in a plastic box
 80 pcs discs and 1 mandren in a plastic box
 240 pcs discs and 4 mandren in a plastic box



RubyPlaton Polishing Strips

Has different kind of strips which has different grin sizes. They are produced combining two different types of grins.

Properties

It allows excellent results in polishing. Flexible and easy to use

- **White** – Extra Fine for Polishing
- **Yellow** – Fine for Finishing
- **Green** – Medium for Countoring
- **Blue** – Coarse for Gross Reduction

Packaging

75 pcs strips in a plastic box
 80 pcs discs and 1 mandren + 75 pcs strips in a plastic box
 240 pcs discs and 4 mandren + 75 pcs strips in a plastic box





RubyPlaton Wooden Wedges

RubyPlaton wooden wedge it has anatomic form for you to make excellent filling.

Properties

- Superthin short
- Thin short
- Thin long
- Medium short
- Medium long

Packaging

100 pcs. wedges in a plastic box
200 pcs. wedges in a plastic box
500 pcs. wedges in a plastic box



RubyPlaton Transparent Wedges

RubyPlaton transparent wedge is formed by special transparent polymer for excellent filling.

Properties

Not easy to break long lived with high quality and durability
Easy to use
Anatomic form
Small size - Medium size

Packaging

100 pcs in a plastic box



Class I



 **WhatsApp**
Get information!



RubyMark Articulating Paper Straight

Great sensitive articulating paper to get simple marking. It is straight form on the 40µ, 80µ and 200µ.

Properties

Non sticky
Smooth surface
Double sided options.

Packaging

40µ – Blue/Red Straight 144 pcs. paper
80µ – Blue/Red Straight 72 pcs. paper
80µ – Blue/Red Straight 144 pcs. paper



Class I



WhatsApp
Get information!



RubyMark Articulating Paper Horseshoe

Great sensitive articulating paper to get simple marking. It is horseshoe form on the 80µ.

Properties

Non sticky
Smooth surface
Double sided options.

Packaging

80µ – Blue/Red Horseshoe 72 pcs. paper

RubyMark Articulating Paper 200 Micron Straight

Great sensitive articulating paper with special fleece paper to get simple marking. It is straight form on the 200µ.

Properties

Non sticky
Smooth surface
Most Accurate Occlusion

Packaging

200µ – Blue Straight 300 pcs. paper
200µ – Red Straight 300 pcs. paper
200µ – Blue Straight 150 pcs. paper – Red Straight 150 pcs. paper



Class I



WhatsApp
Get information!

SOFT START



Rubydent Curing Light

Rubydent Curing Light offers effective and quality curing performance with ergonomic design and metal head.

2 mm resin curing in 1 second
Long life battery / Over 1300 circles under 10 sec mode
Stainless steel head
Replaceable and 360 degree rotating head
Ergonomic design
Black or white colour options

Input	: 100-240VAC 50Hz/60Hz Power
Output	: 5V DC/1A Wave
length	: 385nm-515nm
Light intensity	: 1000mW/cm2– 2500mW/cm2
Battery	: 2600mA

CE Class I



WhatsApp
Get information!



RubyCrownRemover

It is practical crown remover for easy procedure. It removes crown using the way of impact method with arc.

Properties

Can be sterilized at 134 °C
3 types of remover tips
Different types of impact



RubyPlaton Cotton Roll

High quality cotton rolls made from %100 cotton. They are non-linting, non-sterile, soft texture rolls.

Properties

It retracts cheek effective
It does not adhere to mucous membrane
Highly absorbent
Size: 4 x 1 cm
Maintain shape when wet

Packaging

No: 1 / 500 pieces bag
No: 2 / 500 pieces bag
No: 3 / 500 pieces bag

CE Class I




WhatsApp
Get information!


CompacZahn

Acrylic Teeth



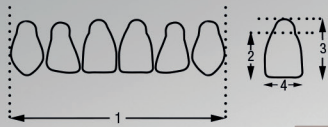

	E1	(1) 43,3 (2) 8,7 (3) 10,4 (4) 8,1
	E2	(1) 42,7 (2) 9,2 (3) 11,0 (4) 8,1
	E3	(1) 45,0 (2) 10,1 (3) 12,1 (4) 8,4
	E4	(1) 48,8 (2) 10,5 (3) 12,1 (4) 8,4
	F2	(1) 42,1 (2) 8,9 (3) 11,6 (4) 7,7
	F3	(1) 42,2 (2) 9,8 (3) 12,4 (4) 7,9
	F4	(1) 48,8 (2) 10,4 (3) 12,5 (4) 8,6
	G1	(1) 40,8 (2) 8,0 (3) 9,7 (4) 7,3
	G2	(1) 44,5 (2) 9,2 (3) 10,5 (4) 8,2
	G3	(1) 44,8 (2) 9,8 (3) 11,6 (4) 8,2
	G4	(1) 48,3 (2) 9,3 (3) 11,8 (4) 8,7
	H1	(1) 40,6 (2) 8,9 (3) 10,4 (4) 7,7
	H2	(1) 42,0 (2) 9,1 (3) 11,4 (4) 7,5
	H3	(1) 42,6 (2) 8,8 (3) 10,4 (4) 7,7
	H4	(1) 43,8 (2) 9,0 (3) 11,4 (4) 8,0

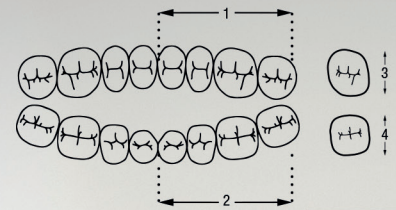

	K1	(1) 32,3 (2) 8,1 (3) 9,7 (4) 4,4
	K3	(1) 33,8 (2) 8,5 (3) 9,7 (4) 4,6
	K5	(1) 34,0 (2) 9,0 (3) 10,9 (4) 5,1
	K6	(1) 35,0 (2) 9,8 (3) 11,0 (4) 5,3
	K8	(1) 40,0 (2) 8,2 (3) 10,1 (4) 6,1
	K9	(1) 38,5 (2) 9,5 (3) 11,2 (4) 6,5
	K10	(1) 40,0 (2) 12,2 (3) 14,0 (4) 5,7
	M1	(1) 29,6 (2) 30,9 (3) 8,5 (4) 7,9
	M2	(1) 29,7 (2) 31,4 (3) 8,6 (4) 8,3
	M3	(1) 31,2 (2) 33,6 (3) 8,9 (4) 8,8
	M4	(1) 31,3 (2) 33,7 (3) 9,2 (4) 8,9
	M5	(1) 33,5 (2) 35,8 (3) 9,8 (4) 9,5

Single Ended



#11000

Double Ended



#11001

Dental Mirror



#11002

Dental Tweezer



#11003

CE Class I



[WhatsApp](#)
Get information!

Pocket Measuring Probe



#11004

Amalgam Carrier



#11005

Matrix Retainer



#11006

Matrix Retainer



#11007

CE Class I



WhatsApp
Get information!

Amalgam Plugger



#11008

Root Canal Plugger

size:1



#11009

Root Canal Plugger

size:2



#11010

Root Canal Plugger

size:3



#11011

CE Class I



WhatsApp
Get information!

Root Canal Spreader

size:1



#11012

Root Canal Spreader

size:2



#11013

Root Canal Spreader

size:3



#11014

Rubber Dam Clamps
for Incisors



#11015

CE Class I



WhatsApp
Get information!

Rubber Dam Clamps
for premolar




#11016

Rubber Dam Clamps
for forceps



#11017

Rubber Dam
frame metal

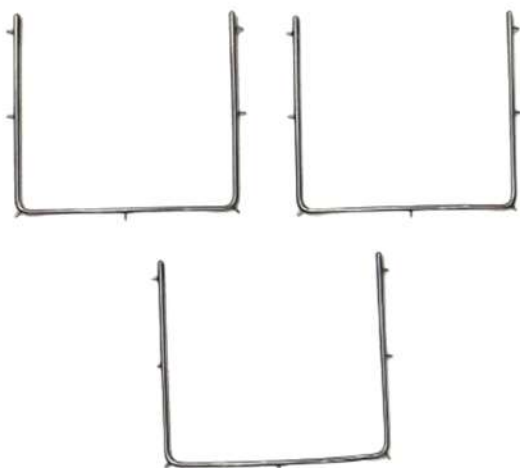

#11018

Endo Clean Stand
endo box



#11019

CE Class I



WhatsApp
Get information!

Endometer Finger Ruler
ring



#11020

Endo Measuring Block



#110021

Rubber Dam Punch



#11022

CE Class I



[WhatsApp](#)
Get information!


ENDOART[®]

FINAL TOUCH



incidental

GLASS SLAB

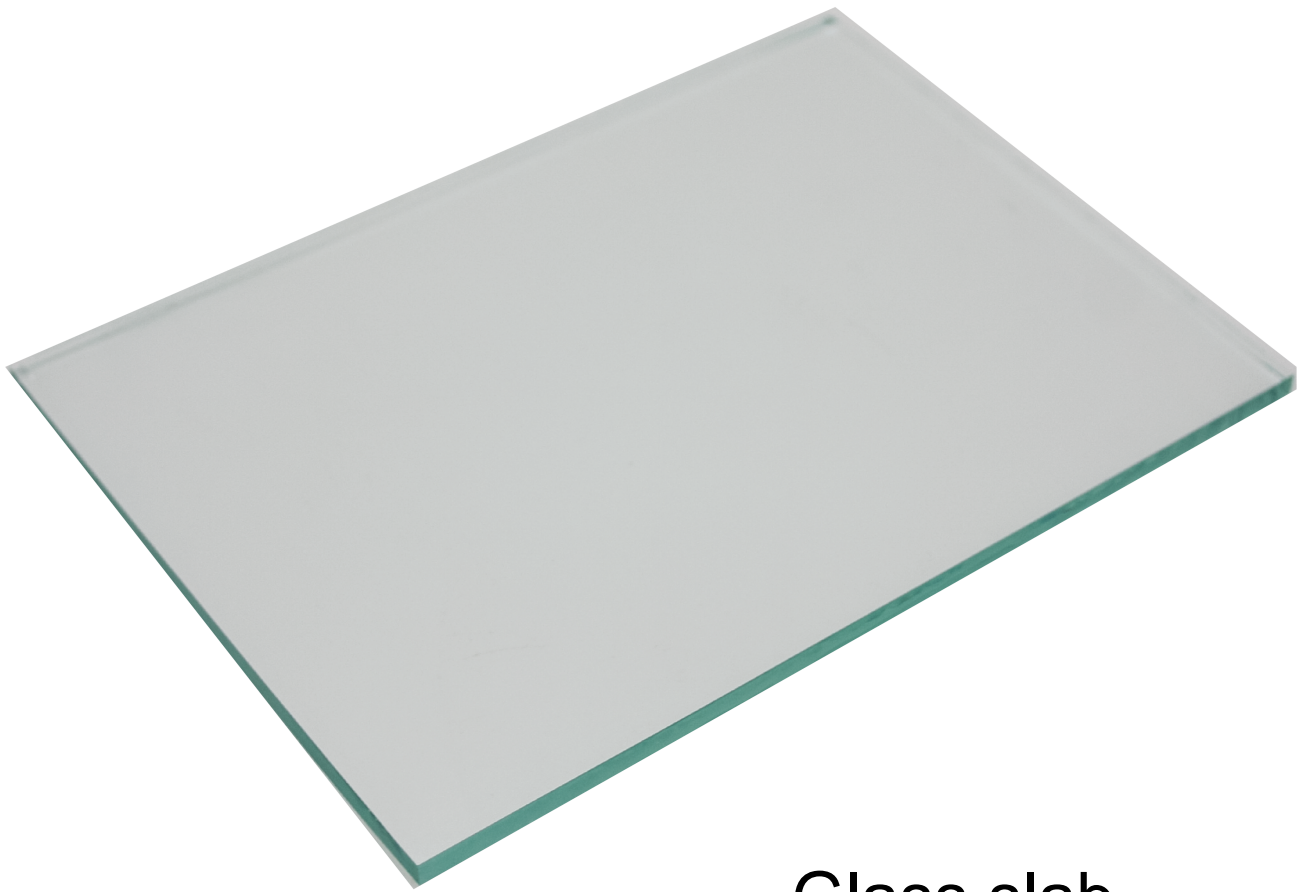

Glass slab
Kod:#11023
Size; 17cmx12cm
Thinck; 6mm

For more information;
export@incidental.com.tr



+905302252254



incidental

İnci Dental Tıbbi Malz. San. ve Tic. Ltd. Şti.

Defterdar Mah. Otakçılar Cad.

E Blok Sit. (Has On Ticaret Merkezi) No:74/4 Eyüpsultan / İstanbul / Turkey

Tel : +90 212 632 63 45 Fax : +90 212 632 97 56

export@incidental.com.tr www.incidental.com.tr